

CANINE COMPLEMENT FACTOR 3

Immunoperoxidase Assay for Determination of C3 in Dog Samples

DIRECTIONS FOR USE

Version2 L6.0 - 4

For Research Use Only, NOT for Diagnostic Purposes

Please Read this Package Insert Completely Before Using This Product

INTENDED USE

The Complement factor 3 (C3) test kits are a highly sensitive two-site enzyme linked immunoassay (ELISA) for measuring C3 in biological samples of dogs.

INTRODUCTION

A number of serum proteins participate in acute inflammatory reactions. These include the complement, coagulation and kinin systems as well as a number of other proteins, known as acute phase proteins, that regulate acute inflammation.

The complement system is a complex set of up to 20 serum proteins. The most abundant and pivotal of the complement components is C3, which has a molecular weight of about 187 kD and consists of an alpha and beta chain.

PRINCIPLE OF THE ASSAY

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the C3 present in samples reacts with the anti-C3 antibodies which have been adsorbed to the surface of polystyrene microtitre wells. After the removal of unbound proteins by washing, anti-C3 antibodies conjugated with horseradish peroxidase (HRP), are added. These enzyme-labeled antibodies form complexes with the previously bound C3. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme varies directly with the concentration of C3 in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of C3 in the test sample. The quantity of C3 in the test sample can be interpolated from the standard curve constructed from the standards, and corrected for sample dilution.

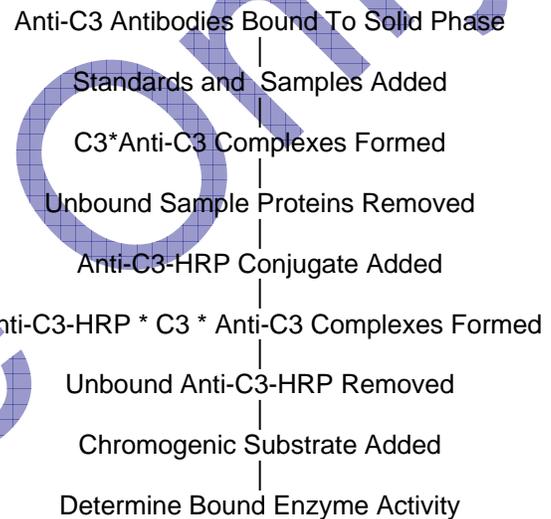


Figure 1.

REAGENTS

(Quantities sufficient for 96 determinations)

1. DILUENT CONCENTRATE (Running Buffer)
One bottle containing 50 ml of a 5X concentrated diluent running buffer.
2. WASH SOLUTION CONCENTRATE
One bottle containing 50 ml of a 20X concentrated wash solution.
3. ENZYME-ANTIBODY CONJUGATE 100X
One vial containing 150 µL of affinity purified anti-Dog C3 antibody conjugated with horseradish peroxidase in a stabilizing buffer.
4. CHROMOGEN-SUBSTRATE SOLUTION
One vial containing 12 mL of 3,3',5,5'-tetramethylbenzidine (TMB) and hydrogen peroxide in citric acid buffer at pH 3.3.

5. STOP SOLUTION

One vial containing 12 ml 0.3 M sulfuric acid.

WARNING: Avoid contact with skin.

6. ANTI-DOG C3 ELISA MICRO PLATE

Twelve removable eight (8) well micro well strips in well holder frame. Each well is coated with affinity purified anti-Dog C3.

7. DOG C3 CALIBRATOR

One vial containing a lyophilized Dog C3 calibrator.

FOR IN VITRO USE ONLY

REAGENT PREPARATION

1. DILUENT CONCENTRATE

The Diluent Solution supplied is a 5X Concentrate and must be diluted 1/5 (1 part buffer concentrate, 4 parts dH₂O) with distilled or deionized water.

2. WASH SOLUTION CONCENTRATE

The Wash Solution supplied is a 20X Concentrate and must be diluted 1/20 with distilled or deionized water (1 part buffer concentrate, 19 part dH₂O). Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.

3. ENZYME-ANTIBODY CONJUGATE

Calculate the required amount of working conjugate solution for each microtitre plate test strip by adding 10 µL Enzyme-Antibody Conjugate to 990 µL of 1X Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.

4. CHROMOGEN-SUBSTRATE SOLUTION

Ready to use as supplied.

5. STOP SOLUTION

Ready to use as supplied.

6. ANTI-DOG C3 ELISA MICRO PLATE

Ready to use as supplied. Unseal Microtiter Pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal. Take clear plastic tape and cover tops of strips to avoid particulates from contaminating wells along with desiccant.

7. DOG C3 CALIBRATOR

Add 1.0 ml of distilled or de-ionized water to the Dog C3 calibrator and mix gently until dissolved. The calibrator is now at a concentration of 10.82 µg/ml

(the reconstituted calibrator should be aliquoted and frozen if future use is intended). Dog C3 standards need to be prepared immediately prior to use (see chart below). Mix well between each step. Avoid foaming.

Standard	ng/ml	Volume added to 1x Diluent	Volume of 1x Diluent
6	200	18.5 µl Dog C3 Calibrator	982 µl
5	100	250 µl standard 6	250 µl
4	50	250 µl standard 5	250 µl
3	25	250 µl standard 4	250 µl
2	12.5	250 µl standard 3	250 µl
1	6.25	250 µl standard 2	250 µl
0	0		500 µl

STORAGE AND STABILITY

The expiration date for the package is stated on the box label.

1. DILUENT

The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4-8°C.

2. WASH SOLUTION

The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions can be stored at room temperature (16-25°C) or at 4-8°C.

3. ENZYME-ANTIBODY CONJUGATE

Undiluted horseradish peroxidase anti-C3 conjugate should be stored at 4-8°C and **diluted immediately prior to use**. The working conjugate solution is stable for up to 8 hours.

4. CHROMOGEN-SUBSTRATE SOLUTION

The Substrate Solution should be stored at 4-8°C and is stable until the expiration date.

5. STOP SOLUTION

The Stop Solution should be stored at 4-8°C and is stable until the expiration date.

6. ANTI-DOG C3 ELISA MICRO PLATE

Anti-Dog C3 coated wells are stable until the expiration date, and should be stored at 4-8°C in sealed foil pouch with desiccant pack.

7. DOG C3 CALIBRATOR

The lyophilized Dog C3 calibrator should be stored at 4°C or frozen until reconstituted. The reconstituted calibrator should be aliquoted out and stored frozen (Avoid multiple freeze-thaw cycles). The working standard solutions should be prepared immediately prior to use and are stable for up to 8 hours.

SPECIMEN COLLECTION AND HANDLING

Blood should be collected by venipuncture and the serum separated from the cells, after clot formation, by centrifugation. Care should be taken to minimize hemolysis, excessive hemolysis can impact your results. Specimens may be shipped at room temperature and then stored refrigerated at 2-8°C if testing is to take place within one week after collection. If testing is to take place later than one week, specimens should be stored at -20°C. Avoid repeated freeze-thaw cycles.

1. Precautions

For any sample that might contain pathogens, care must be taken to prevent contact with open wounds.

2. Additives and Preservatives

No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.

3. Known interfering substances

Azide and thimerosal at concentrations higher than 0.1% inhibits the enzyme reaction.

MATERIAL PROVIDED

See "REAGENTS"

MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipette (2 µL to 200 µL) for making and dispensing dilutions
- Test tubes
- Microtitre washer/aspirator
- Distilled or Deionized H₂O
- Microtitre Plate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer
- Vortex mixer

ASSAY PROTOCOL

DILUTION OF SERUM SAMPLES

The assay for quantification of C3 in serum requires that each test sample be diluted before use. For a single step determination a dilution of serum at 1/40,000 is appropriate for most samples. **For absolute quantification, samples that yield results outside the range of the standard curve, a lesser or greater dilution might be required. If unsure of sample level**

1. To prepare a 1/40,000 dilution of sample, transfer 5 µL of sample to 995 µL of 1X diluent. This gives you a 1/200 dilution. Next, dilute the 1/200 samples by transferring 5 µL, to 995 µL of 1X diluent. You now have a 1/40,000 dilution of your sample. Mix thoroughly at each stage.

PROCEDURE

Bring all reagents to room temperature before use.

1. Add 100 µL of Diluent to each of the wells in 1A & 2A. These will serve for an evaluation of the background associated with the assay.

2. Pipette 100 µL of
Standard 0 (0.0 ng/ml) in duplicate
Standard 1 (6.25 ng/ml) in duplicate
Standard 2 (12.5 ng/ml) in duplicate
Standard 3 (25 ng/ml) in duplicate
Standard 4 (50 ng/ml) in duplicate
Standard 5 (100 ng/ml) in duplicate
Standard 6 (200ng/ml) in duplicate

3. Pipette 100 µL of sample (in duplicate) into pre designated wells.

4. Incubate the micro titer plate at room temperature for twenty (20 ± 2) minutes. Keep plate covered and level during incubation.

5. Following incubation, aspirate the contents of the wells.

6. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with wash buffer, invert the plate then pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual buffer. Repeat 3 times for a total of four washes.

7. Pipette 100 μL of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at room temperature for twenty (20 ± 2) minutes. Keep plate covered and level during incubation.
8. Wash and blot the wells as described in Steps 5/6.
9. Pipette 100 μL of TMB Substrate Solution into each well.
10. Incubate in the dark at room temperature for precisely ten (10) minutes.
11. After ten minutes, add 100 μL of Stop Solution to each well.
12. Determine the absorbance (450 nm) of the contents of each well. Calibrate the plate reader to air.

STABILITY OF THE FINAL REACTION MIXTURE

The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

RESULTS

1. Subtract the average background value from the test values for each sample.

2. Using the results observed for the standards construct a Standard Curve. The appropriate curve fit is that of a four-parameter logistics curve. A second order polynomial (quadratic) or other curve fits may also be used.
3. Interpolate test sample values from standard curve. Correct for sera dilution factor to arrive at the C3 concentration in original sample.

LIMITATION OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.
2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or deionized water, and accuracy of reagent and sample pipettings.

Manufactured by:

Immunology Consultants Laboratory, Inc.
141 N. Elliott Rd | Newberg, OR 97132 | USA
Phone (503) 538-5869 | Fax (503) 554-8510



This document contains information that is proprietary to Immunology Consultants Laboratory. The original recipient of this document may duplicate this document in whole or in part for internal business purposes only, provided that this entire notice appears in all copies. In duplicating any part of this document, the recipient agrees to make every reasonable effort to prevent the unauthorized use and distribution of the proprietary information.