

# MOUSE KIDNEY INJURY MOLECULE (KIM-1)

## Immunoperoxidase Assay for Determination of Kim-1 in Mouse Samples

### DIRECTIONS FOR USE

Version3 1.0 -- 1

For Research Use Only, NOT for Diagnostic Purposes

**Please Read this Package Insert Completely Before Using This Product**

### INTENDED USE

The Kidney Injury Molecule (KIM-1) test kit is a highly sensitive two-site enzyme linked immunoassay (ELISA) for measuring KIM-1 in biological samples of mice.

### INTRODUCTION

KIM-1 is a trans-membrane structural glycoprotein in renal epithelial cells. These cells will shed KIM-1 antigen during regeneration into the urine which makes urinary KIM-1 a useful biomarker for kidney injury.

### PRINCIPLE OF THE ASSAY

The principle of the double antibody sandwich ELISA is represented in Figure 1. In this assay the KIM-1 present in samples reacts with the anti-KIM-1 antibodies which have been adsorbed to the surface of polystyrene microtitre wells. After the removal of unbound proteins by washing, anti-KIM-1 antibodies conjugated with horseradish peroxidase (HRP), are added. These enzyme-labeled antibodies form complexes with the previously bound KIM-1. Following another washing step, the enzyme bound to the immunosorbent is assayed by the addition of a chromogenic substrate, 3,3',5,5'-tetramethylbenzidine (TMB). The quantity of bound enzyme varies directly with the concentration of KIM-1 in the sample tested; thus, the absorbance, at 450 nm, is a measure of the concentration of KIM-1 in the test sample. The quantity of KIM-1 in the test sample can be interpolated from the standard curve constructed from the standards, and corrected for sample dilution.

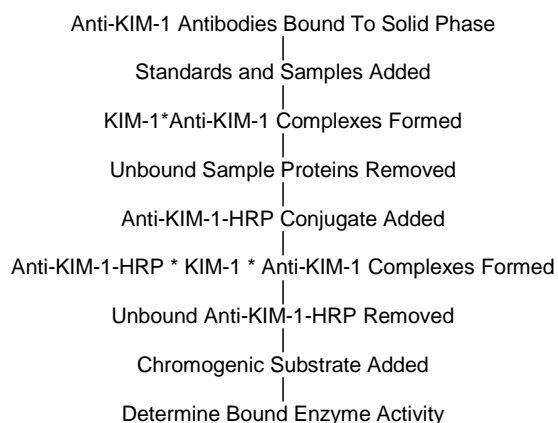


Figure 1.

### REAGENTS (Quantities sufficient for 96 determinations)

- 1. DILUENT (Running Buffer)**  
One bottle containing 50 ml of a 5X concentrated diluent running buffer.
- 2. WASH SOLUTION CONCENTRATE**  
One bottle containing 50 ml of a 20X concentrated wash solution.
- 3. ENZYME-ANTIBODY CONJUGATE 100X**  
One vial containing 150  $\mu$ L of affinity purified anti-Mouse KIM-1 antibody conjugated with horseradish peroxidase in a stabilizing buffer.
- 4. CHROMOGEN-SUBSTRATE SOLUTION**  
One vial containing 12 mL of 3,3',5,5'-tetramethylbenzidine (TMB) and hydrogen peroxide in citric acid buffer at pH 3.3.
- 5. STOP SOLUTION**  
One vial containing 12 ml 0.3 M sulfuric acid.

**WARNING: Avoid contact with skin.**

6. ANTI-MOUSE KIM-1 ELISA MICRO PLATE  
Twelve removable eight (8) well micro well strips in well holder frame. Each well is coated with affinity purified anti-Mouse KIM-1.

7. MOUSE KIM-1 CALIBRATOR  
One vial containing Mouse KIM-1 calibrator.

**FOR IN VITRO USE ONLY**

**REAGENT PREPARATION**

1. DILUENT CONCENTRATE  
The Diluent Solution supplied is a 5X Concentrate and must be diluted 1/5 with distilled or deionized water (1 part buffer concentrate, 4 parts dH<sub>2</sub>O).

2. WASH SOLUTION CONCENTRATE  
The Wash Solution supplied is a 20X Concentrate and must be diluted 1/20 with distilled or deionized water (1 part buffer concentrate, 19 parts dH<sub>2</sub>O). Crystal formation in the concentrate is not uncommon when storage temperatures are low. Warming of the concentrate to 30-35°C before dilution can dissolve crystals.

3. ENZYME-ANTIBODY CONJUGATE  
Calculate the required amount of working conjugate solution for each microtitre plate test strip by adding 10 µL Enzyme-Antibody Conjugate to 990 µL of 1X Diluent for each test strip to be used for testing. Mix uniformly, but gently. Avoid foaming.

4. CHROMOGEN-SUBSTRATE SOLUTION  
Ready to use as supplied.

5. STOP SOLUTION  
Ready to use as supplied.

6. ANTI-MOUSE KIM-1 ELISA MICRO PLATE  
Ready to use as supplied. Unseal Microtiter Pouch and remove plate from pouch. Remove all strips and wells that will not be used in the assay and place back in pouch and re-seal along with desiccant.

7. MOUSE KIM-1 CALIBRATOR  
The Mouse KIM-1 Calibrator should be aliquoted out and stored frozen. The calibrator is now at a concentration of 123.2 ng/ml. **Mouse KIM-1 standards need to be prepared immediately prior to use (see chart below).** Mix well between each step. Avoid foaming.

Standard	ng/ml	Volume added to 1x Diluent	Volume of 1x Diluent
7	5	25 µl Mouse KIM-1 Calibrator	591µl
6	2.5	250 µl standard 7	250 µl
5	1.25	250 µl standard 6	250 µl
4	0.625	250 µl standard 5	250 µl
3	0.3125	250 µl standard 4	250 µl
2	0.156	250 µl standard 3	250 µl
1	0.078	250 µl standard 2	250 µl
0	0		500 µl

**STORAGE AND STABILITY**

The expiration date for the package is stated on the box label.

1. DILUENT  
The 5X Diluent Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions should be stored at 4-8°C.

2. WASH SOLUTION  
The 20X Wash Solution Concentrate is stable until the expiration date. The 1X working solution is stable for at least one week from the date of preparation. Both solutions can be stored at room temperature (16-25°C) or at 4-8°C.

3. ENZYME-ANTIBODY CONJUGATE  
Undiluted horseradish peroxidase anti-KIM-1 conjugate should be stored at 4-8°C and **diluted immediately prior to use.** The working conjugate solution is stable for up to 1 hour when stored in the dark.

4. CHROMOGEN-SUBSTRATE SOLUTION  
The Substrate Solution should be stored at 4-8°C and is stable until the expiration date.

5. STOP SOLUTION  
The Stop Solution should be stored at 4-8°C and is stable until the expiration date.

6. ANTI-MOUSE KIM-1 ELISA MICRO PLATE  
Anti-Mouse KIM-1 coated wells are stable until the expiration date, and should be stored at 4-8°C in sealed foil pouch with desiccant pack.

7. MOUSE KIM-1 CALIBRATOR  
Long Term Storage: Upon receipt, aliquot the calibrator and store them frozen. They will be stable until expiration date. Short Term Storage: the calibrator is stable for up to 10 days at 4°C. The working standard solutions should be prepared immediately prior to use.

## INDICATIONS OF INSTABILITY

If the test is performing correctly, the results observed with the standard solutions should be within 20 % of the expected values.

## SPECIMEN COLLECTION AND HANDLING

Blood should be collected by venipuncture. The serum should be separated from the cells after clot formation by centrifugation. For plasma samples, blood should be collected into a container with an anticoagulant and then centrifuged. Care should be taken to minimize hemolysis, excessive hemolysis can impact your results. For urine samples, collect using standard techniques. The samples should be centrifuged to remove any possible debris. Assay immediately or aliquot and store samples at -20°C. Avoid repeated freeze-thaw cycles.

### 1. Precautions

For any sample that might contain pathogens, care must be taken to prevent contact with open wounds.

### 2. Additives and Preservatives

No additives or preservatives are necessary to maintain the integrity of the specimen. Avoid azide contamination.

### 3. Known interfering substances

Azide and thimerosal at concentrations higher than 0.1% inhibits the enzyme reaction.

## MATERIAL PROVIDED

See "REAGENTS"

## MATERIALS REQUIRED BUT NOT PROVIDED

- Precision pipette (2 µL to 200 µL) for making and dispensing dilutions
- Test tubes
- Microtitre washer/aspirator
- Distilled or Deionized H<sub>2</sub>O
- Microtitre Plate reader
- Assorted glassware for the preparation of reagents and buffer solutions
- Timer

## ASSAY PROTOCOL

### DILUTION OF SAMPLES

The assay for quantification of KIM-1 in samples requires that each test sample be diluted before use. For a single step determination a dilution of 1/3 is

appropriate for most serum, plasma and urine samples. For absolute quantification, samples that yield results outside the range of the standard curve, a lesser or greater dilution might be required. **If unsure of sample level, a serial dilution with one or two representative samples before running the entire plate is highly recommended.**

1. To prepare a 1/3 dilution of sample, transfer 80 µL of sample into 160 µL of 1X diluent. This gives you a 1/3 dilution. Mix thoroughly.

## PROCEDURE

1. **Bring all reagents to room temperature before use.**

2. Pipette 100 µL of

Standard 0 (0.0 ng/ml) in duplicate  
Standard 1 (0.078 ng/ml) in duplicate  
Standard 2 (0.156 ng/ml) in duplicate  
Standard 3 (0.3125 ng/ml) in duplicate  
Standard 4 (0.625 ng/ml) in duplicate  
Standard 5 (1.25 ng/ml) in duplicate  
Standard 6 (2.5 ng/ml) in duplicate  
Standard 7 (5 ng/ml) in duplicate

3. Pipette 100 µL of sample (in duplicate) into pre designated wells.

4. Incubate the micro titer plate at room temperature for one hundred twenty (120 ± 2) minutes. Keep plate covered and level during incubation

5. Following incubation, aspirate the contents of the wells.

6. Completely fill each well with appropriately diluted Wash Solution and aspirate. Repeat three times, for a total of four washes. If washing manually: completely fill wells with wash buffer, invert the plate then pour/shake out the contents in a waste container. Follow this by sharply striking the wells on absorbent paper to remove residual buffer. Repeat 3 times for a total of four washes.

7. Pipette 100 µL of appropriately diluted Enzyme-Antibody Conjugate to each well. Incubate at room temperature for thirty (30 ± 2) minutes. Keep plate covered in the dark and level during incubation.

8. Wash and blot the wells as described in Steps 5/6.

9. Pipette 100 µL of TMB Substrate Solution into each well.

10. Incubate in the dark at room temperature for precisely ten (10) minutes.

11. After ten minutes, add 100  $\mu$ L of Stop Solution to each well.

12. Determine the absorbance (450 nm) of the contents of each well. Calibrate the plate reader to manufacturer's specifications.

### STABILITY OF THE FINAL REACTION MIXTURE

The absorbance of the final reaction mixture can be measured up to 2 hours after the addition of the Stop Solution. However, good laboratory practice dictates that the measurement be made as soon as possible.

### RESULTS

1. Subtract the average background value from the test values for each sample.

2. Using the results observed for the standards construct a Standard Curve. The appropriate curve fit is that of a four-parameter logistics curve. A

second order polynomial (quadratic) or other curve fits may also be used.

3. Interpolate test sample values from standard curve. Correct for sera dilution factor to arrive at the KIM-1 concentration in original samples.

### LIMITATION OF THE PROCEDURE

1. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the information contained in the package insert instructions and with adherence to good laboratory practice.

2. Factors that might affect the performance of the assay include proper instrument function, cleanliness of glassware, quality of distilled or deionized water, and accuracy of reagent and sample pipettings, washing technique, incubation time or temperature.

3. Do not mix or substitute reagents with those from other lots or sources.



Manufactured by:

Immunology Consultants Laboratory, Inc.  
15862 SW 72<sup>nd</sup> Ave  
Suite 150  
Portland, OR 97224 USA  
[www.icllab.com](http://www.icllab.com)

This document contains information that is proprietary to Immunology Consultants Laboratory. The original recipient of this document may duplicate this document in whole or in part for internal business purposes only, provided that this entire notice appears in all copies. In duplicating any part of this document, the recipient agrees to make every reasonable effort to prevent the unauthorized use and distribution of the proprietary information.